

Catalog Number: IR99-409


Anti Caspase-8 cleaved antibody

Package size: 25, 100 μ l
 Store at: -20°C
 MW (kDa): 41-43, 55, 18


Overview

Product Name	Caspase 8 cleaved Asp374 antibody
Product Number	IR99-409
Gene Description	caspase 8, apoptosis-related cysteine peptidase
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human
Recommended Applications Dilutions	Western Blot 1:500 - 1:2000 Immunofluorescence 1:200 - 1:500 Immunohistochemistry (Paraffin) 1:200 - 1:300
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.07 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR99-409 anti-Caspase8 antibody WB image

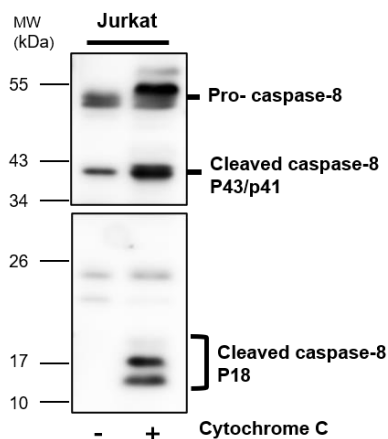

IR99-409 anti-Caspase8 antibody ICC/IF image

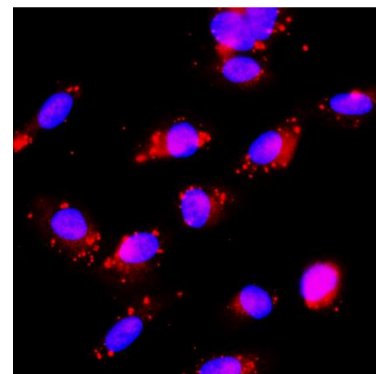


Caspase 8 antibody - at 1/1000 dilution


Lysates at 50 μ g per lane. This blot was produced using a 12% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IR99-409 overnight at 4°C.

Immunofluorescence: cells were fixed with 4% paraformaldehyde for 10 min at RT, permeabilized with 0.1% NP-40 for 10 min at RT then blocked with 5% BSA for 30 min at room temperature. Cells were stained with IR99-409 anti- Caspase 8 antibody (red) at 1:150 and 4°C. DAPI (blue) was used as the nuclear counter stain.

IR99-409 anti-Caspase8 antibody WB image

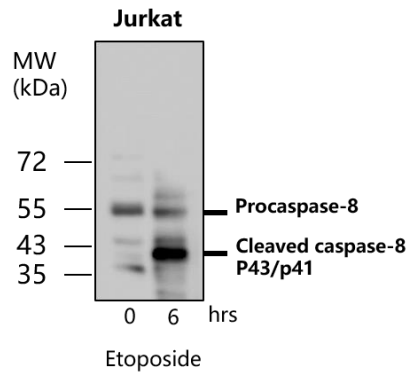

IR99-409 anti-Caspase8 antibody IHC image

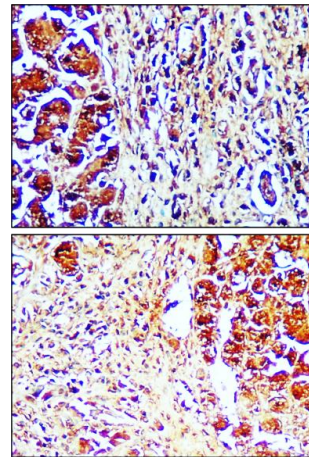


Caspase 8 antibody - at 1/1000 dilution.

Lysates at 50 µg per lane. This blot was produced using a 12% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IR99-409 overnight at 4°C.

Immunohistochemical analysis of paraffin embedded Human cancer tissue labeling Caspase 8 with IR99-409 at 1/200.

Customer feedback Image

IR99-409 anti-Caspase 8 antibody IHC image

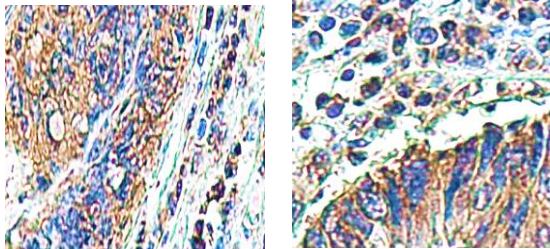

Immunohistochemical of Paraffin-embedded.

Sample: Human colon cancer.
Caspase 8 antibody (IR99-409) : 1-200

IReal Biotechnology Co., Ltd.

 **iReal Antibodies**



**ISO 13485:2016 Quality
Management System
Certified**

 www.irealbio.com

 03-5260005

 irealbio@irealbio.com