

Catalog Number: IR105-435

Anti-Cytochrome C antibody



Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): 14

Overview

Product Name	Cytochrome C antibody
Product Number	IR105-435
Gene Description	cytochrome c, somatic
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, mouse, rat
Recommended Applications Dilutions	Western Blot 1:1000 – 1:2000 Immunofluorescence 1:200 – 1:500 Immunohistochemistry (Paraffin) 1:200 – 1:400
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	0.59 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR105-435 anti-Cytochrome C antibody WB image

IR105-435 anti-Cytochrome C antibody IHC image

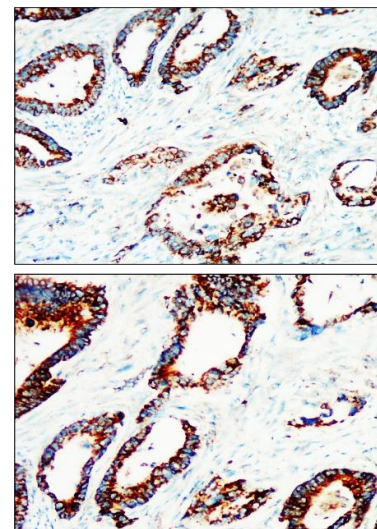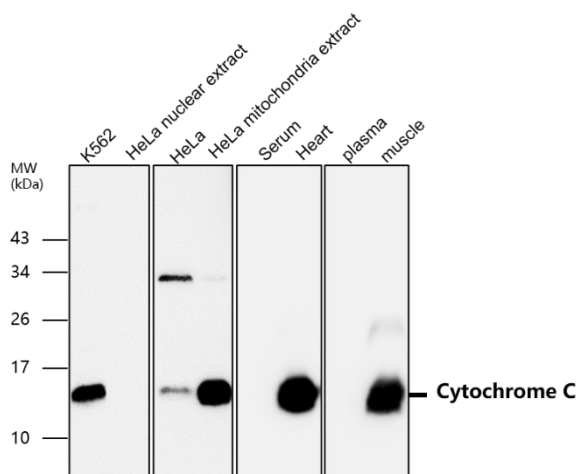


All lanes : Anti- Cytochrome C antibody at 1/1000 dilution
Lysates/proteins at 45 μ g per lane.
This blot was produced using a 15% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IR105-435 overnight at 4°C.

Immunohistochemical analysis of paraffin embedded Human cancer tissue labeling Cytochrome C with IR105-435 at 1/200.

Image

IR105-435 anti-Cytochrome C antibody ICC/IF image

Cytochrome c

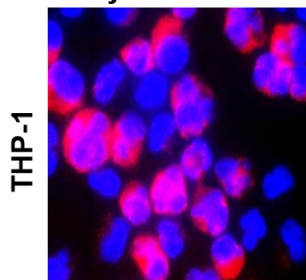


Immunofluorescence: cells were fixed with 4% paraformaldehyde for 10 min at RT, permeabilized with 0.1% NP-40 for 10 min at RT then blocked with 5% BSA for 30 min at room temperature. Cells were stained with IR105-435 anti- Cytochrome C antibody (red) at 1:400 and 4°C. DAPI (blue) was used as the nuclear counter stain.


Customer feedback Image

IR105-435 anti-Cytochrome C antibody IHC image

human colon cancer

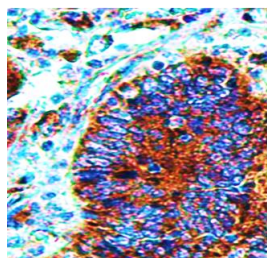

mouse kidney

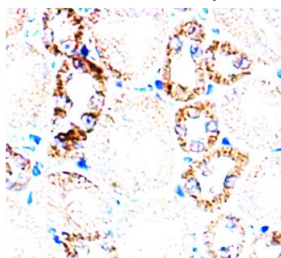

Immunohistochemical analysis of paraffin embedded tissue labeling Cytochrome C with IR105-435 at 1/300.


IReal Biotechnology Co., Ltd.

 **iReal Antibodies**



**ISO 13485:2016 Quality
Management System
Certified**

 www.irealbio.com

 03-5260005

 irealbio@irealbio.com