

Catalog Number: IR307-912


Anti P-selectin / CD62P antibody

Package size: 25, 100 μ l
Store at: -20°C
MW (kDa):


Overview

Product Name	P-selectin / CD62P antibody
Product Number	IR307-912
Gene Description	selectin P (granule membrane protein 140kDa, antigen CD62)
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, mouse, rat
Recommended Applications Dilutions	Western Blot 1:500-1:1000 Immunofluorescence 1:200 – 1:500 Immunohistochemistry (Paraffin) 1:100 – 1:200
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.51 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR307-912 anti P-selectin / CD62P antibody WB image

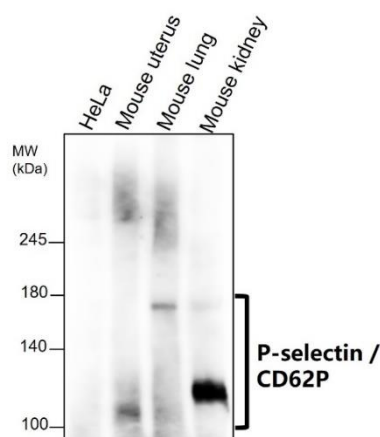


IR307-912 anti P-selectin / CD62P antibody IHC image

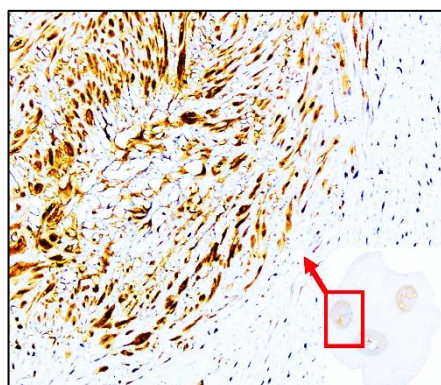

Anti P-selectin / CD62P antibody at 1/500 dilution
Lysates at 50 μ g per lane, This blot was produced using a 5% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IR307-912 overnight at 4°C.


Immunohistochemical analysis of paraffin embedded Human cancer tissue labeling **P-selectin / CD62P** with IR307-912 at 1/100.

IR307-912 anti P-selectin / CD62P antibody IHC image

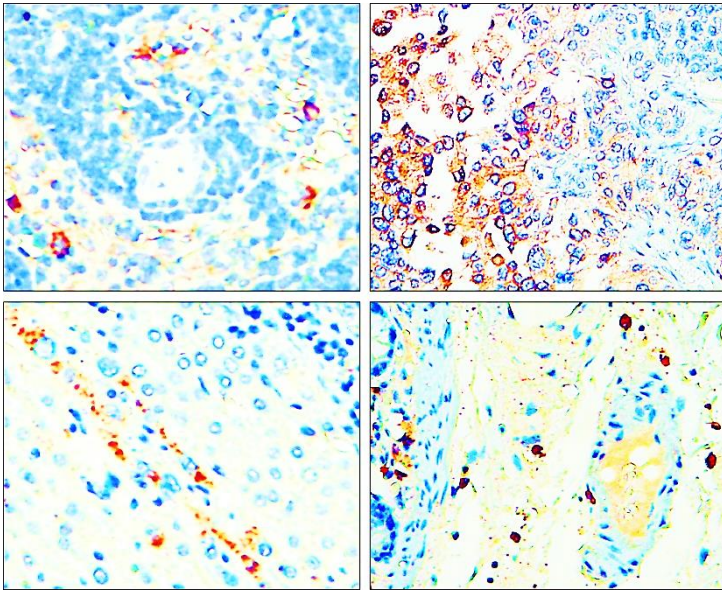

Immunohistochemical analysis of paraffin embedded Human Cancer tissue labeling **P-selectin / CD62P** with IR307-912 at 1/100.

IR307-912 anti P-selectin / CD62P antibody IHC image

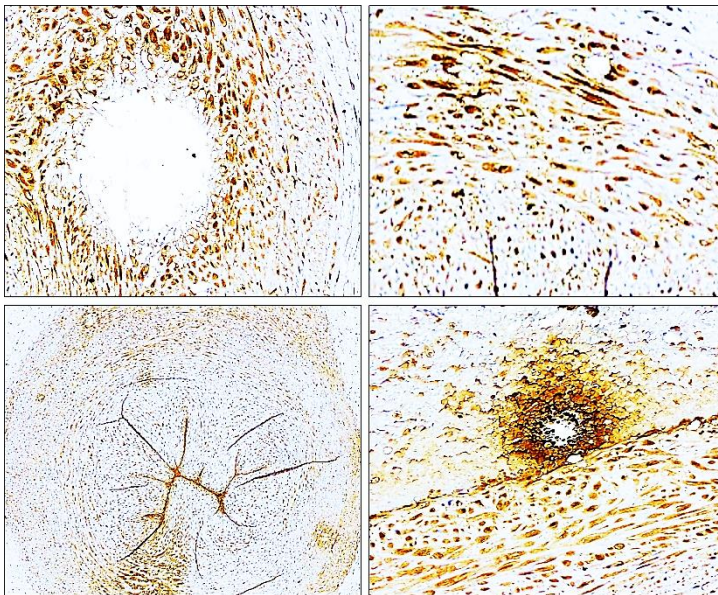

Immunohistochemical analysis of paraffin embedded human tissue labeling **P-selectin / CD62P** with IR307-912 at 1/100 RT 1hr. Epitope Retrieval methods: Citrate Buffer, pH 6.0
[Data from: Li-Tzung Biotechnology. Co](#)

IReal Biotechnology Co., Ltd.

 **iReal Antibodies**



ISO 13485:2016 Quality Management System Certified

www.irealbio.com
[03-5260005](tel:03-5260005)
irealbio@irealbio.com