

Catalog Number: IR306-909


## Anti L-selectin / CD62L antibody

Package size: 25, 100  $\mu$ l  
Store at: -20°C  
MW (kDa): 100-120

### Overview


Product Name	L-selectin / CD62L antibody
Product Number	IR306-909
Gene Description	selectin, lymphocyte
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, mouse, rat
Recommended Applications Dilutions	Western Blot 1:500-1:1000 Immunofluorescence 1:200 – 1:500 Immunohistochemistry (Paraffin) 1:100 – 1:200
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	0.84 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR306-909 anti L-selectin / CD62L antibody WB image

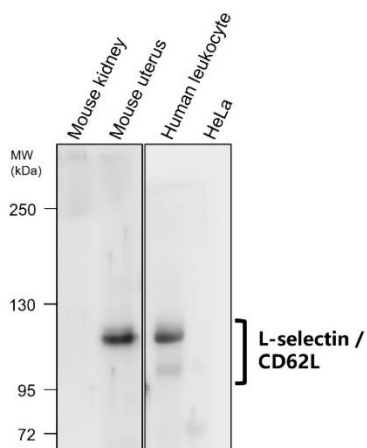


All lanes : Anti L-selectin / CD62L antibody at 1/500 dilution. Lysates at 50  $\mu$ g per lane, This blot was produced using a 5% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IR306-909 overnight at 4°C.

IR306-909 anti L-selectin / CD62L antibody IHC image

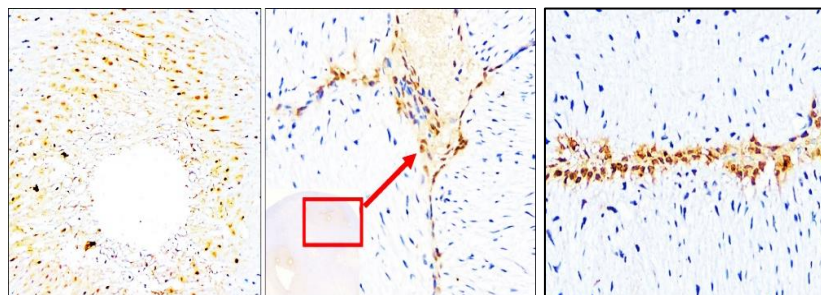

Immunohistochemical analysis of paraffin embedded Human cancer tissue labeling L-selectin / CD62L with IR306-909 at 1/100.

IR306-909 anti L-selectin / CD62L antibody ICC/IF image

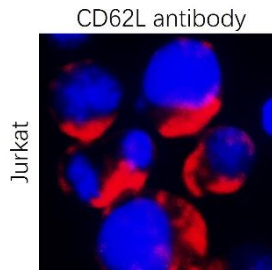

Immunofluorescence: cells were fixed with 4% PFA for 10 min at RT, permeabilized with 0.1% NP-40 for 10 min at RT then blocked with 5% BSA for 30 min at RT. Cells were stained with IR306-909 anti- **L-selectin / CD62** antibody (red) at 1:300 and 4°C. DAPI (blue) was used as the nuclear counter stain.

IReal Biotechnology Co., Ltd.



ISO 13485:2016 Quality  
Management System  
Certified

[www.irealbio.com](http://www.irealbio.com)

03-5260005

[irealbio@irealbio.com](mailto:irealbio@irealbio.com)