

Catalog Number: IR293-964

Anti-Troponin T, cardiac muscle antibody


Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): 39-43

Overview

Product Name	Troponin T, cardiac muscle antibody
Product Number	IR293-964
Gene Description	troponin T, cardiac 3
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, mouse, rat
Recommended Applications Dilutions	Western Blot 1:2000 - 1:4000 Immunohistochemistry (Paraffin) 1:200 – 1:400
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.51 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

Image

IR293-964 anti-Troponin T, cardiac muscle antibody WB image

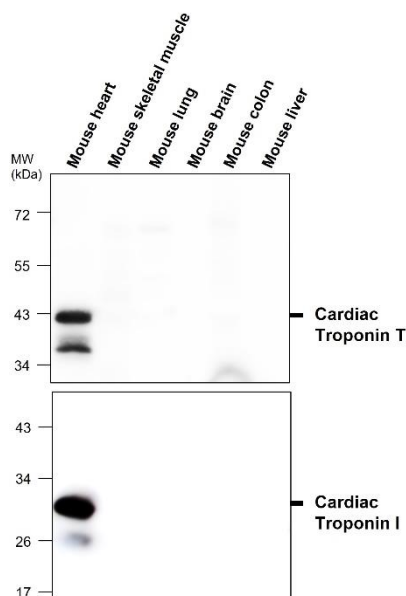

Anti- Troponin T, cardiac muscle antibody - at 1/3000 dilution

Lysates/proteins at 50 μ g per lane
This blot was produced using a 10% SDS-PAGE. The gel was run at 140V for 50 minutes before being transferred onto a Nitrocellulose membrane at 18V for 60 minutes. The membrane was then blocked for an hour before being incubated with IR293-964 overnight at 4°C.

Troponin I, cardiac muscle (IR294-967).