


Catalog Number: IR79-292-PE	Package size: 25, 100 μ l Store at: 4°C
Anti- VCAM-1 / CD106 (PE conjugated) Antibody	

Overview

Product Name	VCAM-1 / CD106 (PE conjugated) antibody
Product Number	IR79-272-PE
Gene Description	vascular cell adhesion molecule 1
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, Mouse
Recommended Applications Dilutions	Immunofluorescence 1:200 – 1:300 Flow Cytometry 1:50 – 1:100
Storage Buffer	PBS (pH 7.4), 0.2% BSA, 0.1% Gelatin, 0.1% NaN3 was added as a preservative
Concentration	1 mg/ml
Purity	Affinity column purified
Storage	This product is stable for 12 months when stored at 4°C. Do not aliquot the antibody. Protect from light. Do not freeze.
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR79-292-PE anti-VCAM-1 / CD106 PE conjugated antibody ICC/IF image

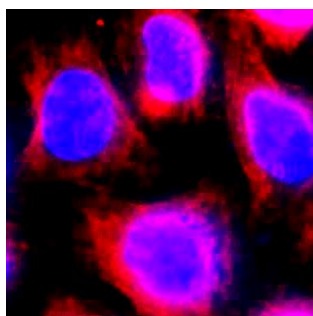

Immunofluorescence: cells were fixed with 4% PFA for 10 min at RT, permeabilized with 0.1% NP-40 for 10 min at RT then blocked with 5% BSA for 30 min at RT. Cells were stained with IR79-292-PE anti-VCAM-1 / CD106 PE conjugated antibody (red) at 1:200 and 4°C. DAPI (blue) was used as the nuclear counter stain.

IReal Biotechnology Co., Ltd.

