

Catalog Number: IRM205


Anti Alpha-1-microglobulin antibody [2H5]

Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): 30-34


Overview

Product Name	Alpha-1-microglobulin antibody [2H5]
Product Number	IRM205
Clonality	Monoclonal
Host	Mouse (mouse IgG1)
Species Reactivity	Human
Recommended Applications Dilutions	Western Blot 1:1000-1:2000 Immunofluorescence 1:200 – 1:500
Storage Buffer	Phosphate-buffered solutions, pH 7.2, 50% Glycerol, 0.025% ProClin 300 was added as a preservative
Concentration	1mg/ml
Purity	Protein A column purified
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IRM205 Alpha-1-microglobulin antibody [2H5] WB image

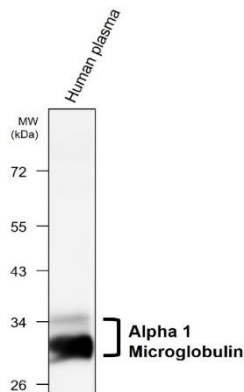

IRM205 Alpha-1-microglobulin antibody [2H5] ICC/IF image

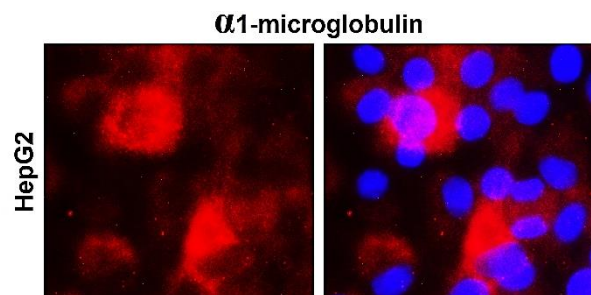

Anti-Alpha-1-microglobulin antibody [2H5] at 1/1000 dilution Lysates at 30 μ g per lane. This blot was produced using a 10% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IRM205 overnight at 4°C.

Immunofluorescence: cells were fixed with 4% paraformaldehyde for 10 min at RT, permeabilized with 0.1% NP-40 for 10 min at RT then blocked with 5% BSA for 30 min at room temperature. Cells were stained with IRM204 anti-**Alpha-1-microglobulin antibody [2H5]** (red) at 1:200 and 4°C. DAPI (blue) was used as the nuclear counter stain.

IReal Biotechnology Co., Ltd.