

Catalog Number: IRST002

Anti-EGFR antibody (Extracellular domain)


Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): 120-180

Overview

Product Name	EGFR antibody (Extracellular domain)
Product Number	IRST002
Gene Description	epidermal growth factor receptor
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human
Recommended Applications Dilutions	Western Blot 1:1000 Immunofluorescence 1:200 – 1:500 Immunohistochemistry (Paraffin) 1:100 – 1:200
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.46 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

Image

IRST002 anti-EGFR antibody WB image

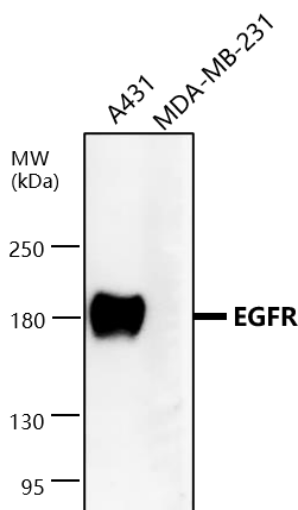

All lanes : Anti-EGFR antibody - at 1/1000 dilution

Lysates/proteins at 50 μ g per lane
This blot was produced using a 5% SDS-PAGE. The gel was run at 140V for 50 minutes before being transferred onto a Nitrocellulose membrane at 18V for 60 minutes. The membrane was then blocked for an hour before being incubated with IRST002 overnight at 4°C.

IReal Biotechnology Co., Ltd.

