

Catalog Number: IR298-1


Anti-COL4A1 / Collagen IV antibody

Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): 200-220


Overview

Product Name	COL4A1 / Collagen IV antibody
Product Number	IR298-1
Gene Description	collagen, type IV
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, mouse, rat
Recommended Applications Dilutions	Western Blot 1:1000 - 1:1500 Immunofluorescence 1:300 – 1:500 Immunohistochemistry (Paraffin) 1:100 – 1:200
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.08 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR298-1 anti-COL4A1 / Collagen IV antibody WB image

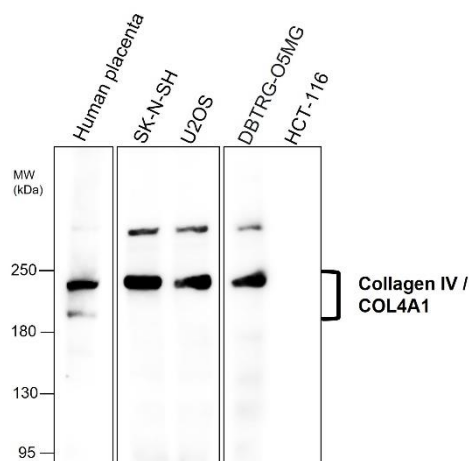


IR298-1 anti-COL4A1 / Collagen IV antibody IHC image

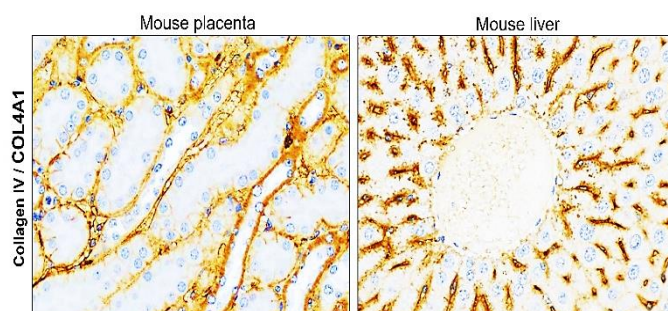

COL4A1 / Collagen IV antibody - at 1/1000 dilution. Lysates at 40 μ g per lane. This blot was produced using a 5% SDS-PAGE. Nitrocellulose membrane was then blocked before being incubated with IR298-1 overnight at 4°C.


Immunohistochemical analysis of paraffin embedded mouse tissue labeling **COL4A1 / Collagen IV** with IR298-1 at 1/100.

IR298-1 anti-COL4A1 / Collagen IV antibody IHC image

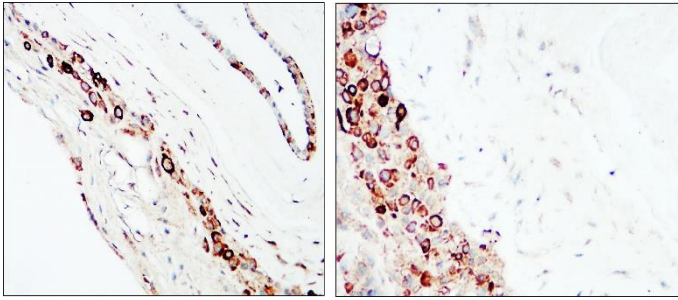

Immunohistochemical analysis of paraffin embedded Human placenta labeling **COL4A1 / Collagen IV** with IR298-1 at 1/100.

IR298-1 anti-COL4A1 / Collagen IV antibody IHC image

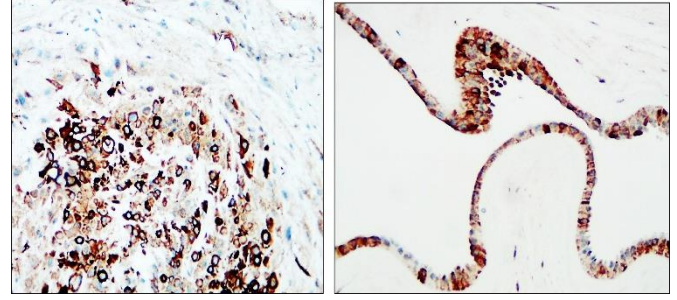

Immunohistochemical analysis of paraffin embedded Human placenta labeling **COL4A1 / Collagen IV** with IR298-1 at 1/100.

IReal Biotechnology Co., Ltd.



ISO 13485:2016 Quality Management System Certified

www.irealbio.com

03-5260005

irealbio@irealbio.com