

Catalog Number: IR304-990


Anti Flotillin-1 antibody

Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): 47
42 (isoform2)


Overview

Product Name	Flotillin-1 antibody
Product Number	IR304-990
Gene Description	flotillin 1
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, mouse, rat
Recommended Applications Dilutions	Western Blot 1:2000 – 1:3000 Immunofluorescence 1:300 – 1:500 Immunohistochemistry (Paraffin) 1:200 – 1:500
Storage Buffer	100mM Tris Glycine, 1% BSA, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	0.96 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR304-990 Flotillin-1 antibody WB image

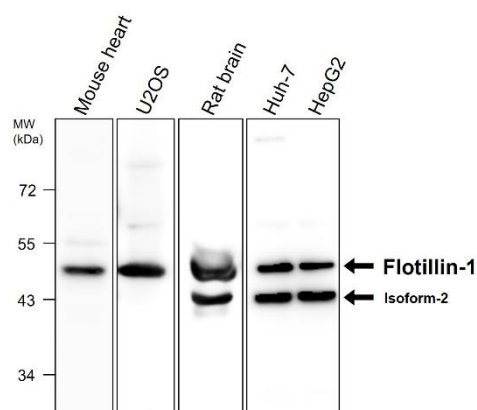


IR304-990 Flotillin-1 antibody ICC/IF image

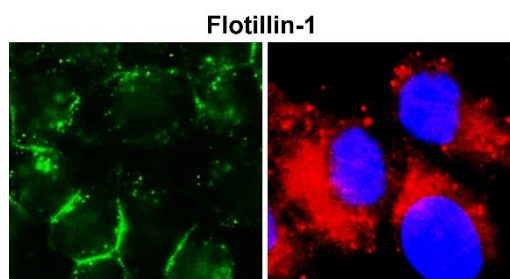

Flotillin-1 at 1/3000 dilution. Lysates at 30 μ g per lane. This blot was produced using a 10% SDS-PAGE. Nitrocellulose membrane was then blocked before being incubated with IR304-990 overnight at 4°C.


Immunohistochemical analysis of paraffin embedded mouse adipose tissue labeling **Flotillin-1** with IR304-990 at 1/500.

IR304-990 Flotillin-1 antibody IHC image

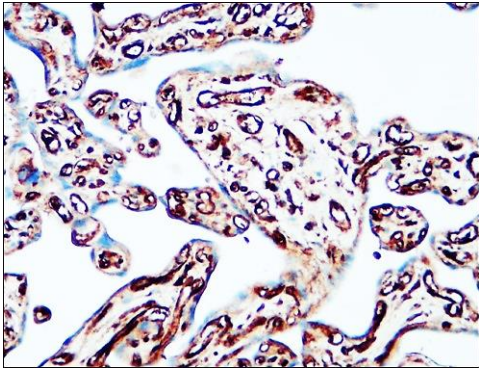

Immunohistochemical analysis of paraffin embedded Human placenta labeling **Flotillin-1** with IR304-990 at 1/100.

IR304-990 Flotillin-1 antibody IHC image

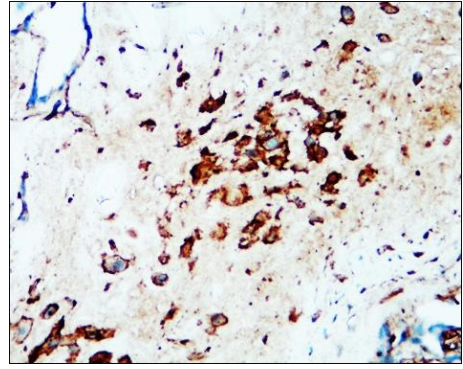

Immunohistochemical analysis of paraffin embedded Human placenta labeling **Flotillin-1** with IR304-990 at 1/100.

IReal Biotechnology Co., Ltd.

 **iReal Antibodies**



**ISO 13485:2016 Quality
Management System
Certified**

 www.irealbio.com

 03-5260005

 irealbio@irealbio.com