

Catalog Number: IRCP220601


Anti-IGLC2 / Immunoglobulin lambda chain antibody

Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): 24-30


Overview

Product Name	IGLC2 / Immunoglobulin lambda chain antibody
Product Number	IR220601
Gene Description	Immunoglobulin lambda constant 2
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human
Recommended Applications Dilutions	Western Blot 1:1000 – 1:2000 Immunofluorescence 1:100 – 1:200
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.5 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IRCP220601 anti-IGLC2 antibody WB image

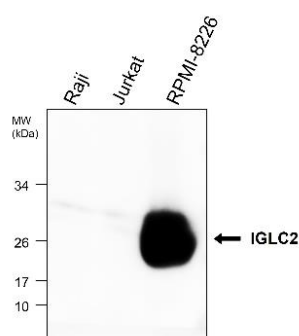

IRCP220601 anti-IGLC2 antibody IHC image

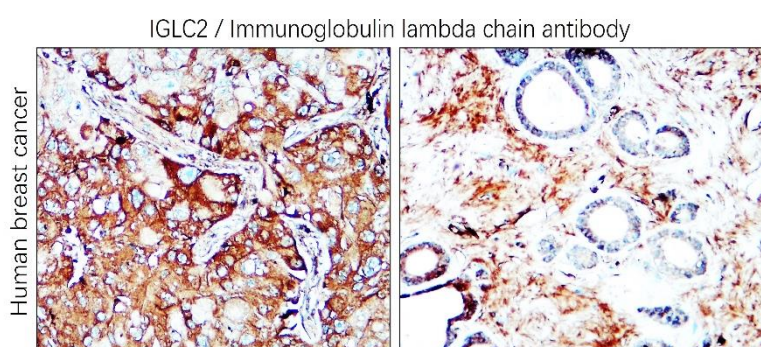


Anti-IGLC2 antibody at 1/1000 dilution
Lysates at 30 μ g per lane. This blot was produced using a 12% SDS-PAGE. Nitrocellulose membrane was then blocked an hour before being incubated with IRCP220601 overnight at 4°C.


Immunohistochemical analysis of paraffin embedded Human breast cancer tissue labeling IGLC2 with IRCP220601 at 1/100.

IRCP220601 anti-IGLC2 antibody IHC image

Human Tonsil/lymphoma

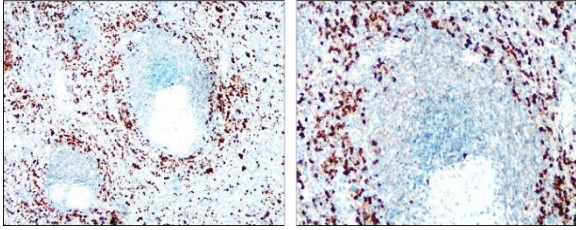


Human Prostate carcinoma

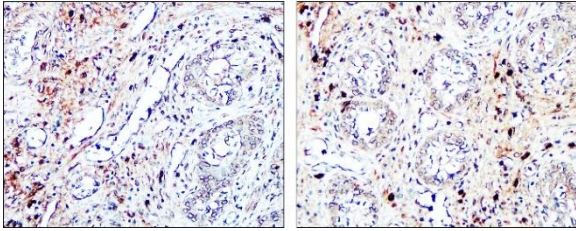

Immunohistochemical analysis of paraffin embedded Human cancer tissue labeling IGLC2 with IRCP220601 at 1/100.

IRCP220601 anti-IGLC2 antibody IHC image

IGLC2 / Immunoglobulin lambda chain antibody

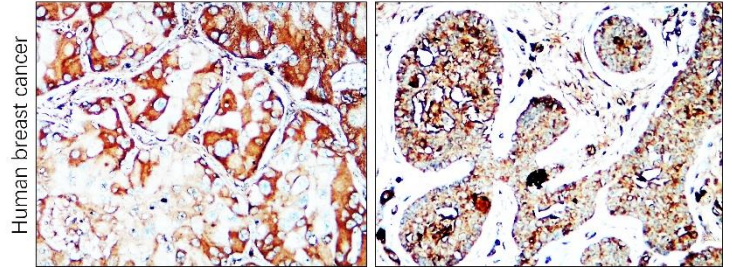

Immunohistochemical analysis of paraffin embedded Human breast cancer tissue labeling IGLC2 with IRCP220601 at 1/100.


IReal Biotechnology Co., Ltd.

 **iReal Antibodies**



**ISO 13485:2016 Quality
Management System
Certified**

 www.irealbio.com

 03-5260005

 irealbio@irealbio.com