

Catalog Number: IR295-2


Anti-Tissue Factor / CD142 antibody

Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): 45-53, >110

Overview

Product Name	Tissue Factor / CD142 antibody
Product Number	IR295-2
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human
Recommended Applications Dilutions	Western Blot 1:1000 - 1:2000 Immunofluorescence 1:300 - 1:600 Immunohistochemistry (Paraffin) 1:100 - 1:300
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.5 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR295-2 anti-Tissue Factor / CD142 antibody WB image

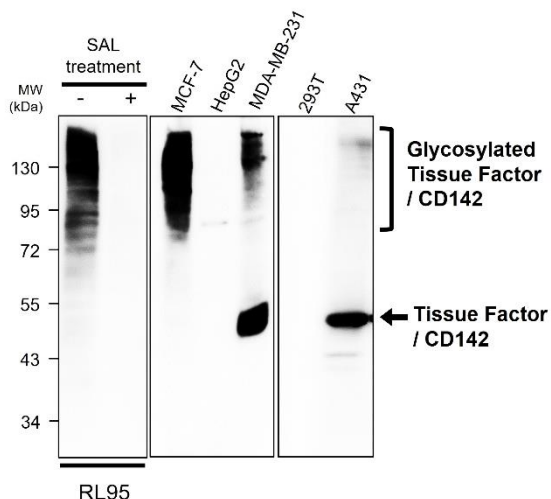

Anti- Tissue Factor / CD142 antibody - at 1/1000 dilution

Lysates/proteins at 30 μ g per lane. This blot was produced using a 10% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IR295-2 overnight at 4°C.