



Catalog No. IRAP5026

Adipogenesis Marker

Antibody Panel

Package: 1 kits (5x25µl)
Store at: -20°C


For Research Use Only. Not For Use In Diagnostic Procedures.

Products Included	Catalog No.	Quantity	Isotype	Application
Adiponectin antibody	IR341-923	25 µl	Rabbit IgG	WB, IF, IHC
Fatty Acid Synthase / FASN Ab	IR342-925	25 µl	Rabbit IgG	WB, IF, IHC
Perilipin-1 antibody	IR343-927	25 µl	Rabbit IgG	WB, IF, IHC
PPAR γ antibody	IR344-930	25 µl	Rabbit IgG	WB, IF, IHC
GAPDH antibody	IR3-8	25 µl	Rabbit IgG	WB, IF, IHC

See www.irealbio.com for individual component applications, species cross-reactivity, dilutions.


Recommended Antibody Dilutions	Western Blot 1:500 - 1:2000 Immunochemistry 1:100 - 1:300 Please visit www.irealbio.com for validation data and a complete listing of recommended companion products.
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Storage / Notes	Store at -20°C. Do not aliquot the antibodies. Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

Adiponectin antibody

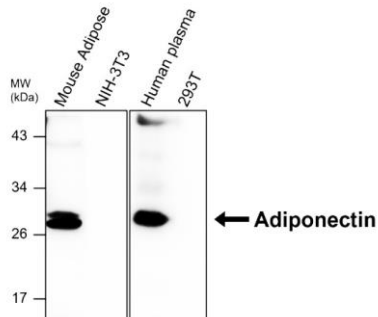


Western blot analysis of cell extracts using **Adiponectin antibody** IR341-923 at 1/500 dilution.

Fatty Acid Synthase / FASN antibody

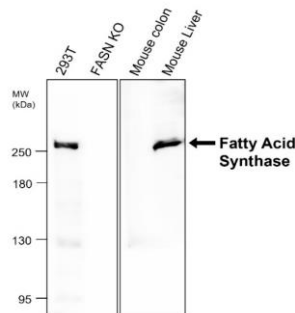


Western blot analysis of cell extracts using **FASN antibody** IR342-925 at 1/1000 dilution.

Perilipin-1 antibody

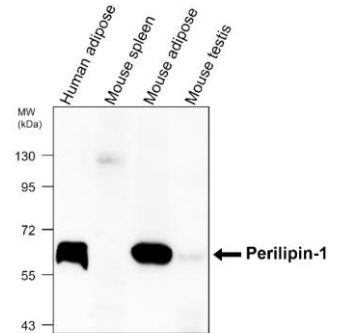


Western blot analysis of cell extracts using **Perilipin-1 antibody** IR343-927 at 1/2000 dilution.

PPAR γ antibody

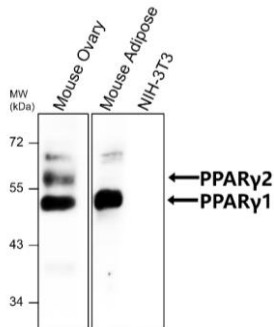

Western blot analysis of cell extracts using **PPAR γ antibody** IR344-930 at 1/1000 dilution.

GAPDH antibody

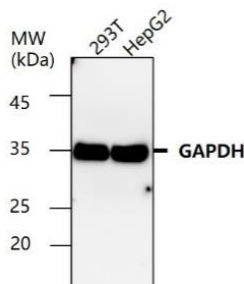

Western blot analysis of cell extracts using **GAPDH antibody** IR3-8 at 1/5000 dilution.

iReal Biotechnology







ISO 13485:2016 Quality Management System Certified



irealbio@irealbio.com

www.irealbio.com

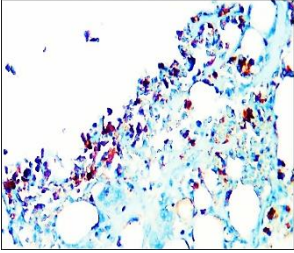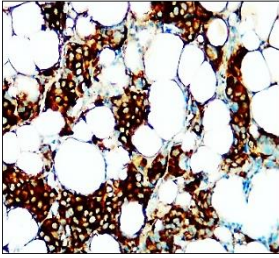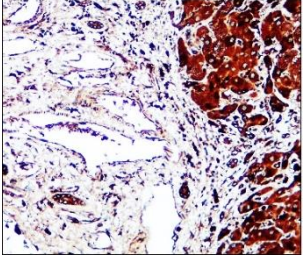

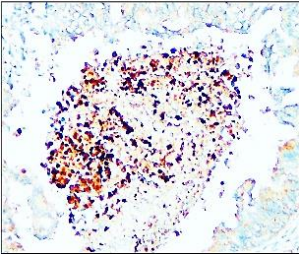Adiponectin antibody	Fatty Acid Synthase / FASN antibody	Perilipin-1 antibody
	<p data-bbox="678 235 869 257">Human breast carcinoma</p> 	<p data-bbox="1244 230 1324 253">Perilipin-1</p> 
<p data-bbox="84 611 488 730">Immunohistochemical analysis of paraffin embedded Human cancer tissue labeling Adiponectin antibody IR341-923 at 1/100.</p>	<p data-bbox="545 627 1023 716">Immunohistochemical analysis of paraffin embedded Human cancer tissue labeling FASN antibody IR342-925 at 1/200.</p>	<p data-bbox="1051 627 1458 716">Immunohistochemical analysis of paraffin embedded Human cancer tissue labeling Perilipin-1 antibody at 1/200.</p>
<p data-bbox="193 766 400 795">PPARγ antibody</p>	<p data-bbox="560 976 951 1021">iReal Biotechnology</p> <p data-bbox="557 1077 888 1171">ISO 13485:2016 Quality Management System Certified</p>  <p data-bbox="564 1196 879 1238">✉ irealbio@irealbio.com</p> <p data-bbox="564 1249 879 1292">🔍 www.irealbio.com</p> 	
<p data-bbox="212 938 397 963">ovarian carcinoma</p> 		
<p data-bbox="84 1339 488 1429">Immunohistochemical analysis of paraffin embedded Human cancer tissue labeling PPARγ antibody IR344-930 at 1/100.</p>		