

WB

IF

IHC

Catalog Number: IRM309

Anti c-Myc antibody [PK-9]


Package size: 25, 100 μ l
Store at: -20°C
MW (kDa):

Overview

Product Name	c-Myc antibody [PK-9]
Product Number	IRM309
Clonality	Monoclonal
Host	Mouse (mouse IgG1)
Species Reactivity	Human
Recommended Applications Dilutions	Western Blot 1:500-1:1000 Immunofluorescence 1:300 – 1:400
Storage Buffer	Phosphate-buffered solutions, pH 7.2, 50% Glycerol, 0.025% ProClin 300 was added as a preservative
Concentration	1mg/ml
Purity	Protein A column purified
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

Image

IRM309 c-Myc antibody [PK-9] WB image

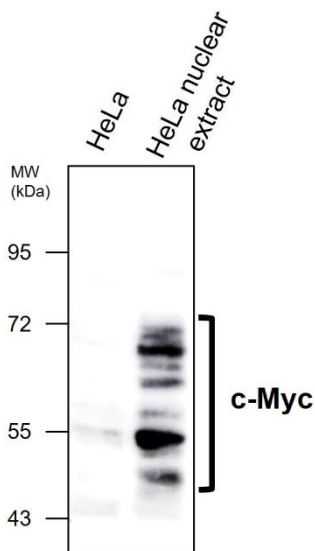

Anti-c-Myc antibody [PK-9] at 1/800 dilution

Lysates at 30 μ g per lane.

This blot was produced using a 10% SDS-PAGE. The gel was run at 140V for 50 minutes before being transferred onto a Nitrocellulose membrane at 18V for 60 minutes. The membrane was then blocked for an hour before being incubated with IRM309 overnight at 4°C.

