

Catalog Number: IRM022


Anti-p23 monoclonal antibody [G-7]

Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): 53

Overview

Product Name	p23 monoclonal antibody [G-7]
Product Number	IRM022
Clonality	Monoclonal
Host	Mouse (mouse IgG1)
Species Reactivity	Human, Mouse, Rat
Recommended Applications Dilutions	Western Blot 1:500-1:1000 Immunofluorescence 1:200 – 1:300
Storage Buffer	Phosphate-buffered solutions, pH 7.2, 50% Glycerol, 0.025% ProClin 300 was added as a preservative
Concentration	0.5mg/ml
Purity	Protein A column purified
Storage / Notes	Long time storage is recommended at -20°C. Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IRM022 p23 monoclonal antibody [G-7] WB image

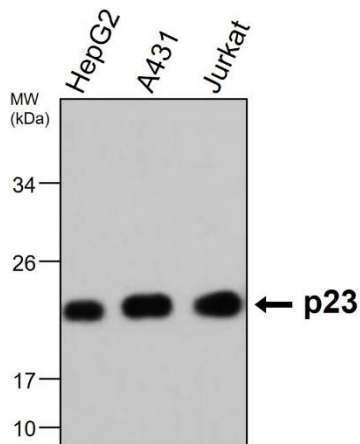

Anti- p23 monoclonal antibody [G-7] at 1/500 dilution. Lysates at 30 μ g per lane. NC membrane was then blocked for an hour before being incubated with IRM022 overnight at 4°C.

