

Catalog Number: IRM303

## Anti-CD21 / CR2 antibody [PC-3]


Package size: 25, 100 µl  
Store at: -20°C  
MW (kDa): 165

### Overview

Product Name	CD21 / CR2 antibody [PC-3]
Product Number	IRM303
Clonality	Monoclonal
Host	Mouse (mouse IgG1)
Species Reactivity	Human
Recommended Applications Dilutions	Western Blot 1:500-1:1000 Immunofluorescence 1:300 – 1:500
Storage Buffer	Phosphate-buffered solutions, pH 7.2, 50% Glycerol, 0.025% ProClin 300 was added as a preservative
Concentration	1mg/ml
Purity	Protein A column purified
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

### Image

#### IRM303 CD21 / CR2 antibody [PC-3] WB image

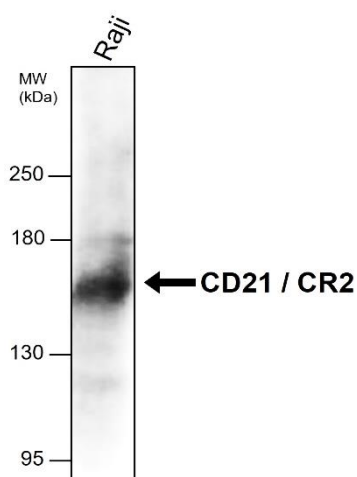

#### Anti-CD21 / CR2 antibody [PC-3] at 1/500 dilution

Lysates at 30 µg per lane.  
This blot was produced using a 5% SDS-PAGE. The gel was run at 140V for 50 minutes before being transferred onto a Nitrocellulose membrane at 18V for 60 minutes. The membrane was then blocked for an hour before being incubated with IRM303 overnight at 4°C.

iReal Biotechnology Co., Ltd.