

Catalog Number: IR361

## Anti D-Dimer antibody


Package size: 25, 100  $\mu$ l  
Store at: -20°C  
MW (kDa): -

### Overview

|                                    |   |
|------------------------------------|---|
| Product Name                       | D-Dimer antibody  |
| Product Number                     | IR361   |
| Clonality                          | Polyclonal  |
| Host                               | Rabbit  |
| Species Reactivity                 | Human   |
| Recommended Applications Dilutions | Western Blot 1:500 – 1:1000<br>Immunofluorescence 1:200 – 1:400   |
| Storage Buffer                     | 100mM Tris Glycine, 30% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative                              |
| Concentration                      | 1.14 mg/ml  |
| Purity                             | Affinity column purified  |
| Storage                            | Store at +4°C for short term storage. Long time storage is recommended at -20°C                                     |
| Notes                              | Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user. |

### IR361 anti D-Dimer antibody WB image

IR361  
D-Dimer antibody

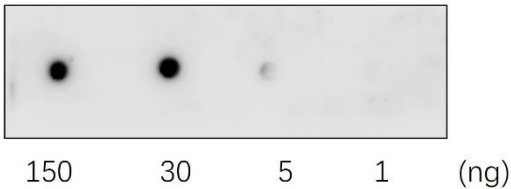

#### Anti D-Dimer antibody at 1/1000 dilution

Specificity and sensitivity of **D-Dimer** Antibody (IR361)

Nitrocellulose membrane was then blocked for an hour before being incubated with IR361 overnight at 4°C.

iReal Biotechnology Co., Ltd.

