

Catalog Number: IRM304

Anti-CD31 / PECAM1 antibody [PD-4]


Package size: 25, 100 μ l
Store at: -20°C
MW (kDa):

Overview

Product Name	CD31 / PECAM1 antibody [PD-4]
Product Number	IRM304
Clonality	Monoclonal
Host	Mouse (mouse IgG1)
Species Reactivity	Human
Recommended Applications Dilutions	Western Blot 1:500 Immunofluorescence 1:200 – 1:300 Immunohistochemistry (Paraffin) 1:100 – 1:200
Storage Buffer	Phosphate-buffered solutions, pH 7.2, 50% Glycerol, 0.025% ProClin 300 was added as a preservative
Concentration	1mg/ml
Purity	Protein A column purified
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IRM304 CD31 / PECAM1 antibody [PD-4] IHC image

Human placenta

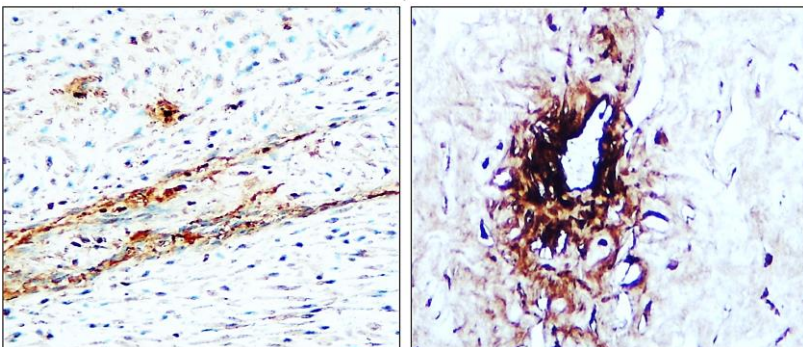

Immunohistochemical analysis of paraffin embedded Human placenta tissue labeling **CD31 / PECAM1** antibody with IRM304 at 1/100.

