

Catalog Number: IR351-953

Anti-CDX2 antibody

Package size: 25, 100 μ l

Store at: -20°C

MW (kDa): 30-38



Overview

Product Name	CDX2 antibody
Product Number	IR351-953
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, Mouse, Rat
Recommended Applications Dilutions	Western Blot 1:1000 - 1:2000 Immunofluorescence 1:300 - 1:600 Immunohistochemistry (Paraffin) 1:100 - 1:200
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.5 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR351-953 anti-CDX2 antibody IHC image

IR351-953 anti-CDX2 antibody WB image

Human colorectal adenocarcinoma

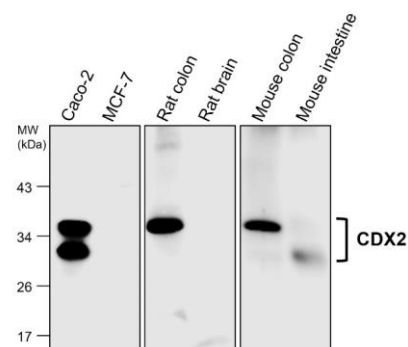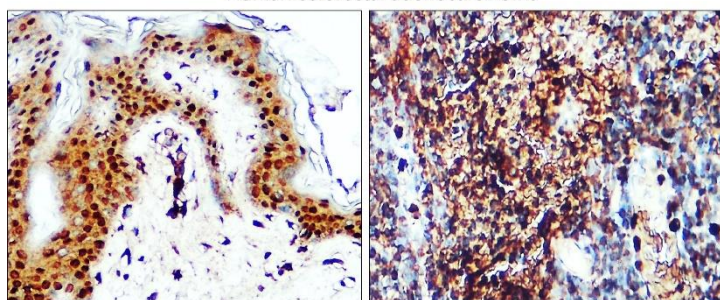

Immunohistochemical analysis of paraffin embedded Human colorectal adenocarcinoma tissue labeling CDX2 with IR351-953 at 1/100.

Anti-CDX2 antibody - at 1/1000 dilution Lysates at 35 μ g per lane. This blot was produced using a 12% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IR351-953 overnight at 4°C.

IReal Biotechnology Co., Ltd.