

Catalog Number: IRM024


Anti-NOX2 / gp91phox, monoclonal antibody

Package size: 25, 100 μ l
 Store at: -20°C
 MW (kDa):


Overview

Product Name	NOX2 / gp91phox, monoclonal antibody
Product Number	IRM024
Clonality	Monoclonal
Host	Mouse (mouse IgG1)
Species Reactivity	Human
Recommended Applications Dilutions	Western Blot 1:500 Immunofluorescence 1:100 – 1:200 Immunohistochemistry (Paraffin) 1:50 – 1:100
Storage Buffer	Phosphate-buffered solution, pH 7.2, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	0.5 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IRM024 NOX2 gp91phox, monoclonal antibody
 WB image

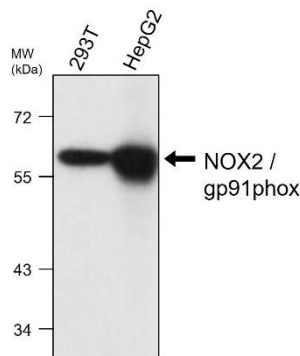

IRM024 NOX2 gp91phox, monoclonal antibody
 IHC image

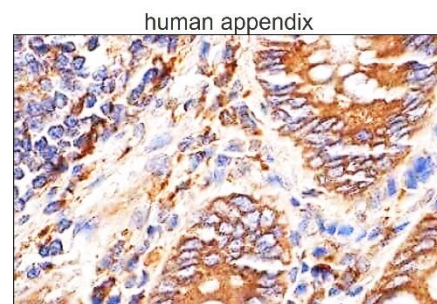

NOX2 gp91phox, monoclonal antibody at 1/500 dilution.
 Lysates/ at 40 μ g per lane. Nitrocellulose membrane was then blocked for an hour before being incubated with IRM024 overnight at 4°C.

Immunoperoxidase staining of formalin fixed, paraffin-embedded human appendix tissue.

