

Catalog Number: IRM028


Anti-Hexokinase II / HK2, monoclonal antibody

Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): 120


Overview

Product Name	Hexokinase II / HK2, monoclonal antibody
Product Number	IRM028
Clonality	Monoclonal
Host	Mouse (mouse IgG1)
Species Reactivity	Human
Recommended Applications Dilutions	Western Blot 1:500 – 1:1000 Immunofluorescence 1:100 – 1:200 Immunohistochemistry (Paraffin) 1:50 – 1:100
Storage Buffer	Phosphate-buffered solution, pH 7.2, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	0.5 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.


IRM028 Hexokinase II / HK2, monoclonal antibody WB image

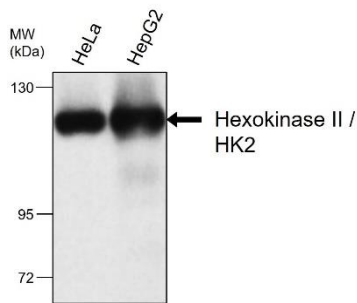

IRM028 Hexokinase II / HK2, monoclonal antibody IHC image

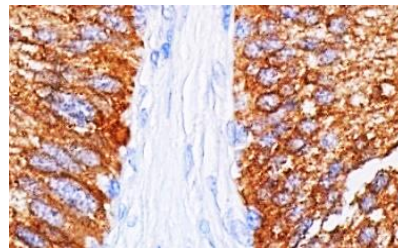

IRM028 Hexokinase II / HK2, monoclonal antibody ICC/IF image

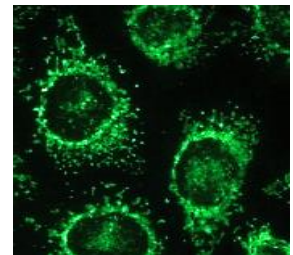

Hexokinase II / HK2 monoclonal antibody at 1/500 dilution. Lysates/ at 35 μ g per lane. Nitrocellulose membrane was then blocked for an hour before being incubated with IRM028 overnight at 4°C.

Immunoperoxidase staining of formalin fixed, paraffin-embedded human epididymis tissue.

Immunofluorescence of cell.
Sample: HeLa
Green: Hexokinase II / HK2 antibody (IRM028): 1-100

IReal Biotechnology Co., Ltd.