

WB
IF
IHC

Catalog Number: IRM031


Anti-STAT1, monoclonal antibody

Package size: 25, 100 μ l
 Store at: -20°C
 MW (kDa): 84, 91


Overview

Product Name	STAT1, monoclonal antibody
Product Number	IRM031
Clonality	Monoclonal
Host	Mouse (mouse IgG1)
Species Reactivity	Human, Mouse, Rat
Recommended Applications Dilutions	Western Blot 1:500 – 1:1000 Immunofluorescence 1:100 – 1:200 Immunohistochemistry (Paraffin) 1:50 – 1:100
Storage Buffer	Phosphate-buffered solution, pH 7.2, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	0.5 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IRM031 STAT1, monoclonal antibody WB image

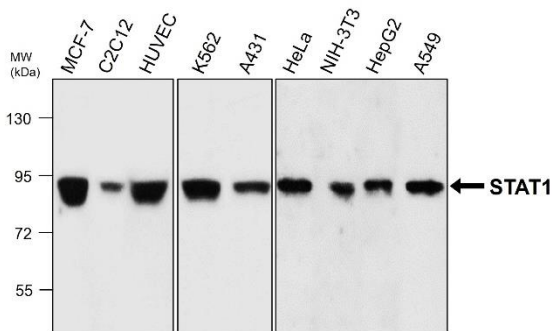

IRM031 STAT1, monoclonal antibody ICC/IF image

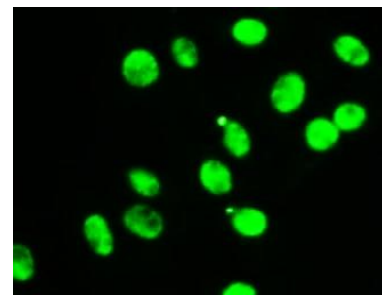

STAT1, monoclonal antibody at 1/500 dilution. Lysates/ at 35 μ g per lane. Nitrocellulose membrane was then blocked for an hour before being incubated with IRM031 overnight at 4°C.

Immunofluorescence of cell.

Sample: HeLa cells, serum-starved overnight and treated with Human Interferon- α 1 (hIFN- α 1) 30 min
 Green: STAT1, monoclonal antibody (IRM031): 1-200

