

Catalog Number: IRCP220803

Anti-ALPK2 antibody

Package size: 25, 100 μ l


Store at: -20°C

MW (kDa): 235


Overview

Product Name	ALPK2 antibody
Product Number	IRCP220803
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, Mouse, Rat
Recommended Applications Dilutions	Western Blot 1:1000 Immunofluorescence 1:200 – 1:500 Immunohistochemistry (Paraffin) 1:200 – 1:500
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.5 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IRCP220803 anti-ALPK2 antibody WB image

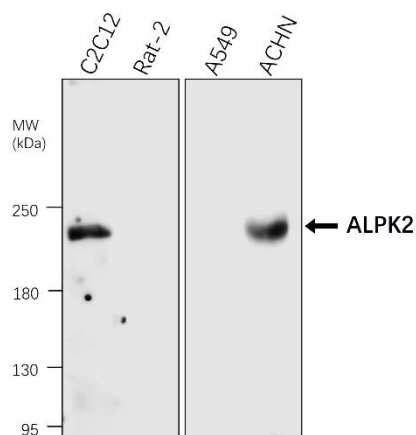

IRCP220803 anti-ALPK2 antibody IHC image

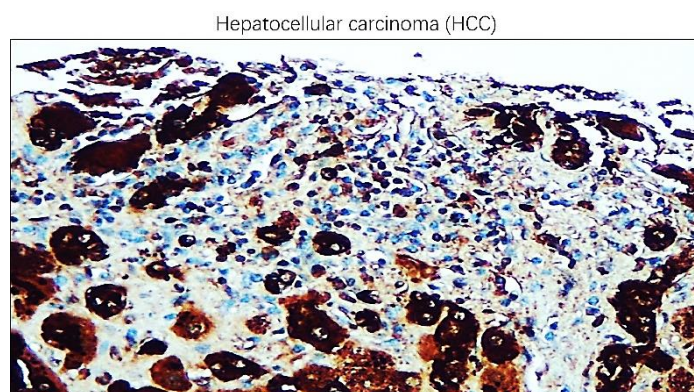


Anti-ALPK2 antibody - at 1/1000 dilution
Lysates/proteins at 40 μ g per lane. This blot was produced using a 5% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IR220803 overnight at 4°C.

Immunohistochemical analysis of paraffin embedded Hepatocellular carcinoma (HCC) tissue labeling ALPK2 antibody with IRCP220803 at 1/200.

IRCP220803 anti-ALPK2 antibody IHC image

Hepatocellular carcinoma (HCC)

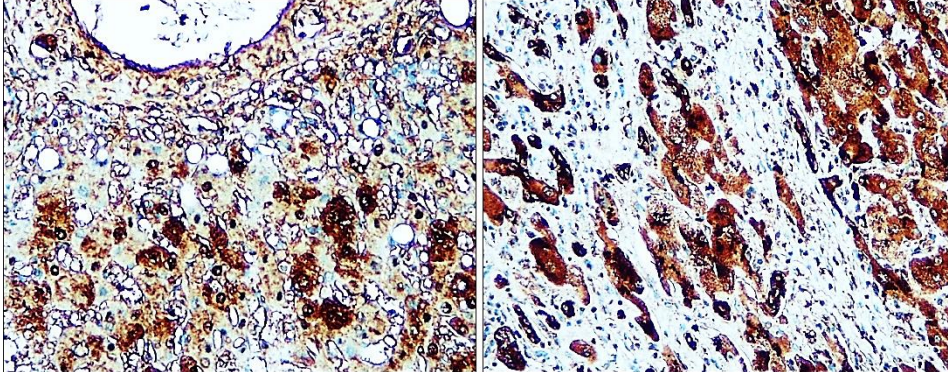

Immunohistochemical analysis of paraffin embedded Hepatocellular carcinoma (HCC) tissue labeling ALPK2 antibody with IRCP220803 at 1/200.

IReal Biotechnology Co., Ltd.

 **iReal Antibodies**



**ISO 13485:2016 Quality
Management System
Certified**

 www.irealbio.com

 03-5260005

 irealbio@irealbio.com