

Catalog Number: IR341-923


Anti-Adiponectin antibody

Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): 27-28


Overview

Product Name	Adiponectin antibody
Product Number	IR341-923
Gene Description	adiponectin, C1Q and collagen domain containing
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, Mouse, Rat
Recommended Applications Dilutions	Western Blot 1:500 – 1:1000 Immunofluorescence 1:200 – 1:300 Immunohistochemistry (Paraffin) 1:100 – 1:200
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.11 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR341-923 anti-Adiponectin antibody WB image

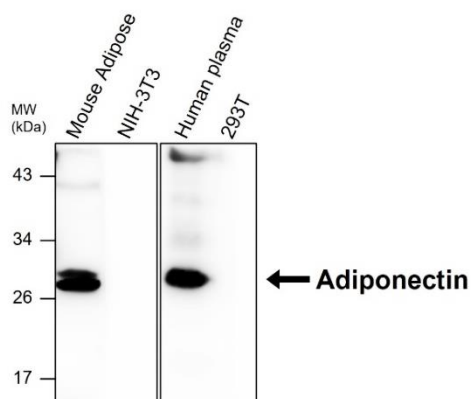

IR341-923 anti-Adiponectin antibody IHC image

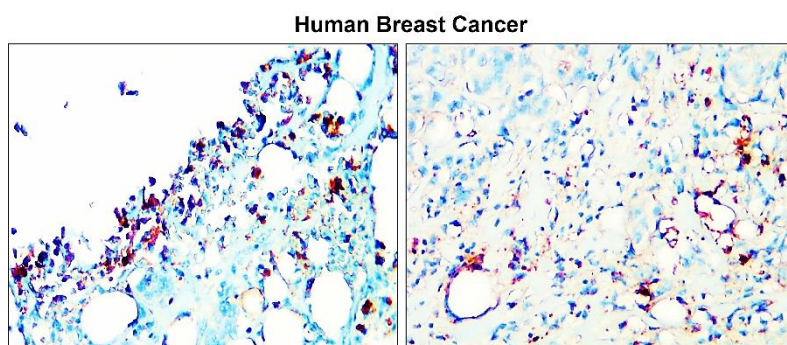

Anti-Adiponectin antibody at 1/1000 dilution
Lysates at 50 μ g per lane. This blot was produced using a 12% SDS-PAGE. Nitrocellulose membrane was then blocked an hour before being incubated with IR341-923 overnight at 4°C.

Immunohistochemical analysis of paraffin embedded Human cancer tissue labeling Adiponectin with IR341-923 at 1/100.

IReal Biotechnology Co., Ltd.