

WB

IP

ELISA

Catalog Number: IRT008

Anti GST-tag [6D] monoclonal antibody

Package size: 25, 100 μ l
Store at: -20°C
MW (kDa):


Overview

Product Name	GST-tag [6D] monoclonal antibody
Product Number	IRT008
Clonality	Monoclonal
Host	Mouse
Isotype	Mouse IgG1
Immunogen	GST recombinant protein
Recommended Applications Dilutions	Western Blot 1:2000-1:5000 ELISA 1:5000 – 1:10000
Storage Buffer	50mM Tris Glycine, 40% Glycerol (pH7.4). 0.025% ProClin 300 was added as a preservative.
Concentration	1 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

Image

IRT008 anti GST-tag [6D] monoclonal antibody WB image

GST recombinant protein

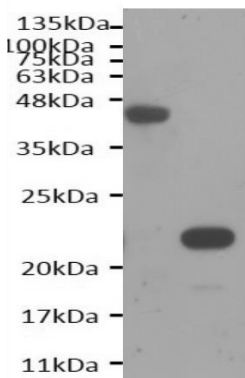

All lanes : Anti GST-tag [6D] monoclonal antibody at 1/3000 dilution

This blot was produced using a 15% SDS-PAGE. The gel was run at 140V for 50 minutes before being transferred onto a Nitrocellulose membrane at 18V for 60 minutes. The membrane was then blocked for an hour before being incubated with IRT008 overnight at 4°C.

iReal Biotechnology Co., Ltd.

