

Catalog Number: IR68-229


## Anti-MMP3 antibody

Package size: 25, 100  $\mu$ l  
 Store at: -20°C  
 MW (kDa): 60

### Overview

Product Name	MMP3 antibody
Product Number	IR68-229
Gene Description	matrix metalloproteinase 3
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human
Recommended Applications Dilutions	Western Blot 1:500 Immunofluorescence 1:100 – 1:200
Storage Buffer	100mM Tris Glycine, 1% BSA, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.5 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

### Image

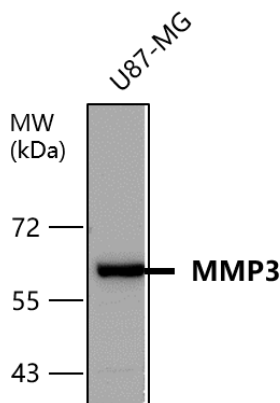

IR68-229 anti-MMP3 antibody WB image

All lanes : Anti-MMP3 antibody at 1/500 dilution

Lysates/proteins at 50  $\mu$ g per lane  
 This blot was produced using a 10% SDS-PAGE. The gel was run at 140V for 50 minutes before being transferred onto a Nitrocellulose membrane at 18V for 60 minutes. The membrane was then blocked with 3% BSA for an hour before being incubated with IR68-229 overnight at 4°C.

IReal Biotechnology Co., Ltd.