



Catalog No. IRAP5035

Mitochondria marker Antibody

Panel

Package: 1 kits (5x25µl)
Store at: -20°C


For Research Use Only. Not For Use In Diagnostic Procedures.

Products Included	Catalog No.	Quantity	Isotype	Application
Cytochrome C antibody	IR105-435	25 µl	Rabbit IgG	WB, IF, IHC
VDAC1 antibody	IR16-2	25 µl	Rabbit IgG	WB, IF, IHC
BAX antibody	IR93-389	25 µl	Rabbit IgG	WB, IF, IHC
AIF monoclonal antibody	IRM405	25 µl	Mouse IgG1	WB, IF, IHC
GAPDH antibody	IR3-8	25 µl	Rabbit IgG	WB, IF, IHC

See www.irealbio.com for individual component applications, species cross-reactivity, dilutions.


Recommended Antibody Dilutions	Western Blot 1:500 - 1:1000 Immunochemistry 1:100 - 1:300 Please visit www.irealbio.com for validation data and a complete listing of recommended companion products.
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Storage / Notes	Store at -20°C. Do not aliquot the antibodies. Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

Cytochrome C antibody

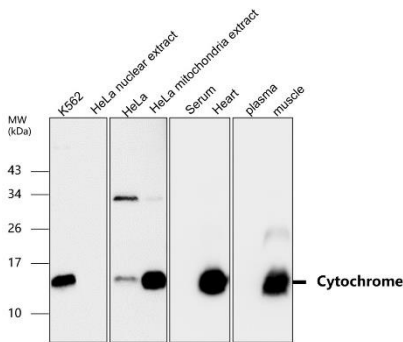


Western blot analysis of cell extracts using Cytochrome C IR105-435 at 1/1000 dilution.

VDAC1 antibody

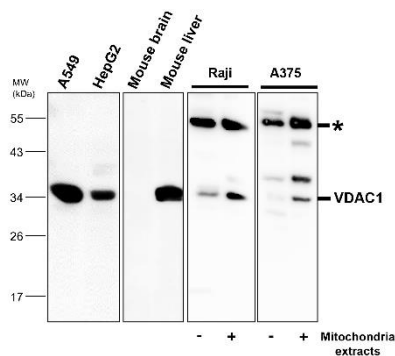


Western blot analysis of cell extracts using VDAC1 IR16-2 at 1/1000 dilution.

BAX antibody

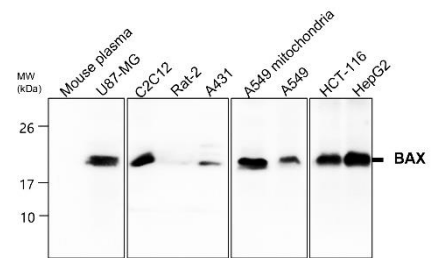


Western blot analysis of cell extracts using BAX IR93-389 at 1/1000 dilution.

AIF monoclonal antibody

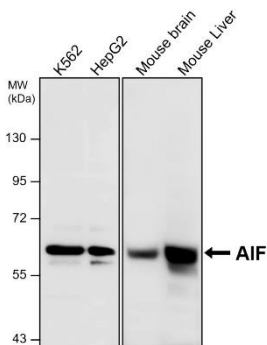

Western blot analysis of cell extracts using AIF monoclonal IRM405 at 1/1000 dilution.

GAPDH antibody

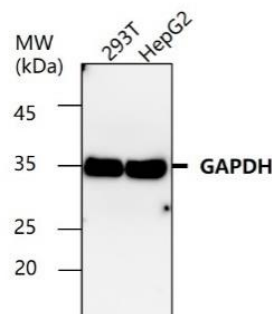

Western blot analysis of cell extracts using GAPDH antibody IR3-8 at 1/5000 dilution.

iReal Biotechnology

ISO 13485:2016 Quality Management System Certified



irealbio@irealbio.com

www.irealbio.com

