

Catalog Number: IRR005


Anti-mouse CD45 antibody [30-F11]

Package size: 25,100 μ l
Store at: -20°C
MW (kDa): -

Overview


Product Name	mouse CD45 antibody [30-F11]
Product Number	IRR005
Clonality	Monoclonal
Host	Rat (Rat IgG2a)
Species Reactivity	Mouse
Recommended Applications Dilutions	Immunohistochemistry 1:100 Flow Cytometry 1:50 – 1:100 Immunofluorescence 1:200 – 1:300
Storage Buffer	Phosphate-buffered solution, pH 7.2, containing 0.025% ProClin 300, 20% Glycerol
Concentration	0.5 mg/ml
Purity	Affinity column purified
Storage	Store at -20°C.
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IRR005 mouse CD45 antibody [30-F11] IHC image

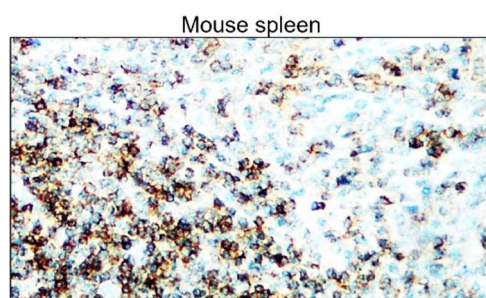


Immunohistochemical analysis of paraffin embedded mouse spleen tissue labeling mouse CD45 antibody [30-F11] with IRR005 at 1/200 RT 1hr.

IRR005 mouse CD45 antibody [30-F11] FACS image

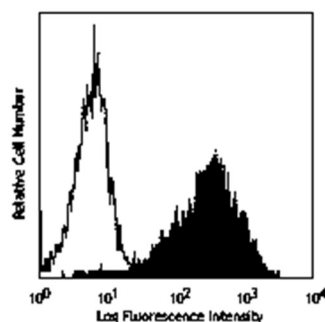

Mouse peripheral blood monocytes intracellularly stained with IRR005 mouse CD45 [30-F11] monoclonal antibody

IRR005 mouse CD45 antibody [30-F11] ICC/IF image

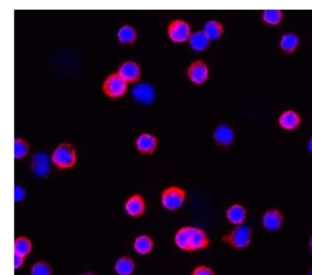

CD45 was detected in immersion fixed mouse splenocytes using Rat Anti-Mouse CD45 Monoclonal Antibody (IRR005) at 1-100, at room temperature. counterstained with DAPI (blue).

iReal Biotechnology Co., Ltd.

