

Catalog Number: IRPA2303

Anti SARS-CoV-2 (COVID-19)


Spike protein S2 antibody

Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): 180-360

Overview

Product Name	SARS-CoV-2 (COVID-19) Spike protein S2 antibody
Product Number	IRPA2303
Clonality	Polyclonal
Host	Rabbit
Species	Virus
Recommended Applications Dilutions	Western Blot 1:2000 – 1:14000 Immunofluorescence 1:500 – 1:1000 Neutralizing Reported in the literature, not verified in house.
Storage Buffer	PBS (Ph7.4), 20% Glycerol. 0.025% ProClin 300 was added as a preservative
Concentration	1 mg/ml
Purity	Protein A column purified
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IRPA2303 anti-SARS-CoV-2 (COVID-19) Spike protein S2 antibody WB image

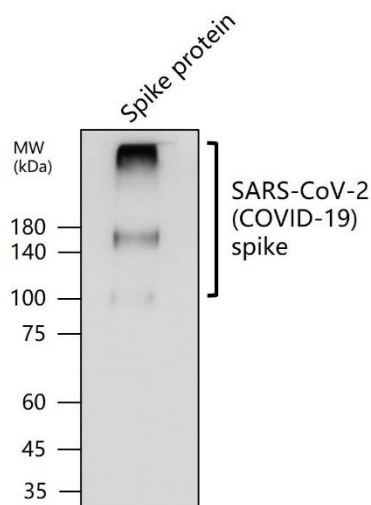

Western blot analysis of **SARS-CoV-2 (COVID-19) Spike protein S2 antibody**.

Proteins at 0.1 μ g per lane.

SARS-CoV-2 (COVID-19) Spike protein was separated on SDS-PAGE gels under reducing conditions. After transferring the proteins to onto NC membranes.