

Catalog Number: IRM042

Anti SARS-CoV-2 (COVID-19)


Spike RBD antibody [10F]

Package size: 25, 100 μ l
Store at: -20°C
MW (kDa):

Overview

Product Name	SARS-CoV-2 (COVID-19) Spike RBD antibody [10F]
Product Number	IRM042
Clonality	Monoclonal
Host	Mouse (mouse IgG1)
Species	Virus
Recommended Applications Dilutions	Western Blot 1:1000 – 1:12000 Immunofluorescence 1:500 – 1:1000 Neutralizing Reported in the literature, not verified in house.
Storage Buffer	PBS (Ph7.4), 20% Glycerol. 0.025% ProClin 300 was added as a preservative
Concentration	1 mg/ml
Purity	Protein A column purified
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IRM042 anti-SARS-CoV-2 (COVID-19) Spike RBD antibody [10F] WB image

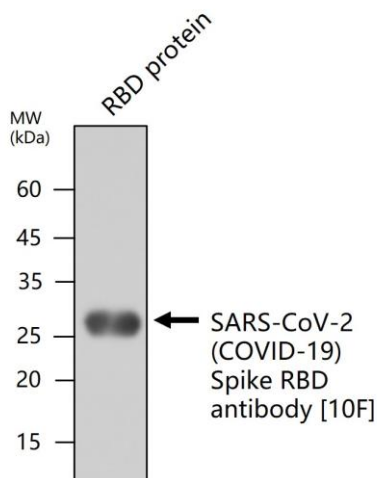

Western blot analysis of **SARS-CoV-2 (COVID-19) Spike RBD antibody [10F]**.

Proteins at 0.1 μ g per lane.

SARS-CoV-2 (COVID-19) Spike RBD protein was separated on SDS-PAGE gels under reducing conditions. After transferring the proteins to onto NC membranes.