

Catalog Number: IRM301

Anti-BDNF antibody [PA-1]


Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): 32-34, 15

Overview

Product Name	BDNF antibody [PA-1]
Product Number	IRM301
Clonality	Monoclonal
Host	Mouse (mouse IgG1)
Species Reactivity	Human
Recommended Applications Dilutions	Western Blot 1:500 – 1:1000 Immunofluorescence 1:300 – 1:500
Storage Buffer	Phosphate-buffered solutions, pH 7.2, 50% Glycerol, 0.025% ProClin 300 was added as a preservative
Concentration	1mg/ml
Purity	Protein A column purified
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

Image

IRM301 BDNF antibody [PA-1] WB image

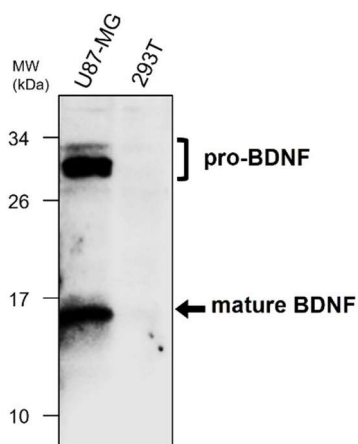

Anti-BDNF antibody [PA-1] at 1/800 dilution

Lysates at 30 μ g per lane.

This blot was produced using a 15% SDS-PAGE. The gel was run at 140V for 50 minutes before being transferred onto a Nitrocellulose membrane at 18V for 60 minutes. The membrane was then blocked for an hour before being incubated with IRM301 overnight at 4°C.

