

Catalog Number: IRM321


## Anti IL1 beta antibody [MD5]

Package size: 25, 100  $\mu$ l  
Store at: -20°C  
MW (kDa): 36


### Overview

Product Name	IL1 beta antibody [MD5]
Product Number	IRM321
Clonality	Monoclonal
Host	Mouse (mouse IgG1)
Species Reactivity	Human
Recommended Applications Dilutions	Western Blot 1:500 Immunofluorescence 1:200 – 1:500 Immunohistochemistry (Paraffin) 1:100 – 1:200
Storage Buffer	Phosphate-buffered solutions, pH 7.2, 50% Glycerol, 0.025% ProClin 300 was added as a preservative
Concentration	1mg/ml
Purity	Protein A column purified
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

### IRM321 IL1 beta antibody [MD5] WB Image

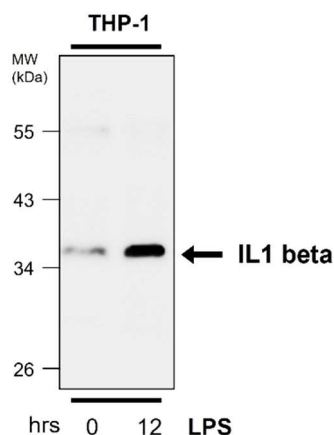

### IRM321 IL1 beta antibody [MD5] IHC Image

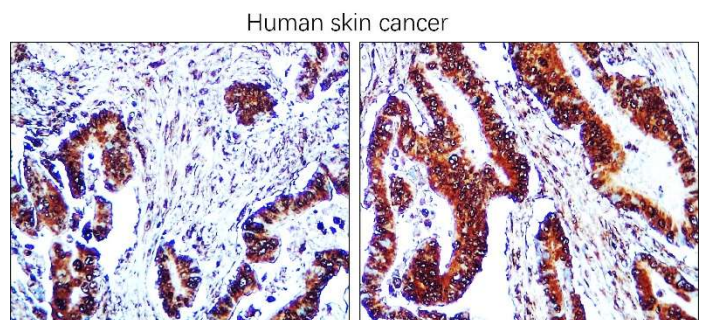

**Anti-IL1 beta antibody [MD5] antibody at 1/500.** Lysates at 35  $\mu$ g per lane. This blot was produced using a 10% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IRM321 overnight at 4°C.

Immunohistochemical analysis of paraffin embedded Human Hepatocellular Carcinoma and placenta tissue labeling **IL1 beta antibody [MD5]** antibody with IRM321 at 1/200.

iReal Biotechnology Co., Ltd.