

Catalog Number: IRFS4042


Mouse anti-rabbit IgG, FITC Conjugated

Package size: 100 µl
Store at: 4°C


Overview

Product Name	Mouse anti-rabbit IgG, FITC Conjugated
Product Number	IRFS4042
Label	FITC Conjugated
Clonality	Monoclonal
Host	Mouse
Species Reactivity	Rabbit IgG
Recommended Applications Dilutions	Immunofluorescence 1:100 – 1:400 Immunohistochemistry 1:100 – 1:400 flow cytometry (0.5-1 µg per 1 x 10 ⁶ cells) Optimal dilution to be determined by titration. (Other Assays User Optimized)
Storage Buffer	PBS containing 1% stabilizer protein and 0.02% sodium azide.
Concentration	0.4 mg/ml
Storage	Store at 4° C, **DO NOT FREEZE**. Stable for one year from the date of shipment. Non-hazardous. No MSDS required.
SOURCE	mouse anti-rabbit IgG-FITC is an affinity purified secondary antibody raised in mouse against rabbit IgG and conjugated to FITC (fluorescein isothiocyanate).

Image

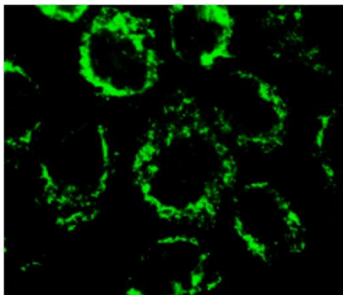

Immunofluorescence staining of formalin-fixed A431 cells showing mitochondrial localization. Detection reagent used: mouse anti-rabbit IgG, FITC Conjugated, IRFS4042.

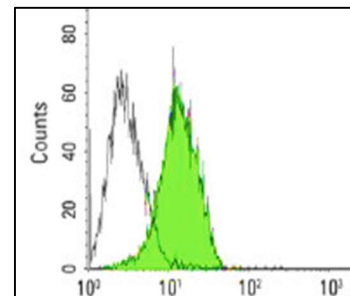

FCM analysis of fixed and permeabilized Raji cells stained with CD45 antibody (IR77-263), followed by FITC-conjugated mouse anti-rabbit IgG: IRFS4042. Black line histogram represents the isotype control.

IReal Biotechnology Co., Ltd.