

Catalog Number: IRFS4052


## Mouse anti-rabbit IgG, PE Conjugated

Package size: 100 µl  
Store at: 4°C


### Overview

Product Name	Mouse anti-rabbit IgG, PE Conjugated
Product Number	IRFS4052
Label	PE Conjugated
Clonality	Monoclonal
Host	Mouse
Species Reactivity	Rabbit IgG
Recommended Applications Dilutions	Immunofluorescence 1:100 – 1:400 Immunohistochemistry 1:100 – 1:400 flow cytometry (0.5-1 µg per 1 x 10 <sup>6</sup> cells) Optimal dilution to be determined by titration. (Other Assays User Optimized)
Storage Buffer	PBS containing 1% stabilizer protein and 0.02% sodium azide.
Concentration	0.4 mg/ml
Storage	Store at 4° C, <b>**DO NOT FREEZE**</b> . Stable for one year from the date of shipment. Non-hazardous. No MSDS required.
SOURCE	mouse anti-rabbit IgG-PE is an affinity purified secondary antibody raised in mouse against rabbit IgG and conjugated to PE (phycoerythrin).

### Image

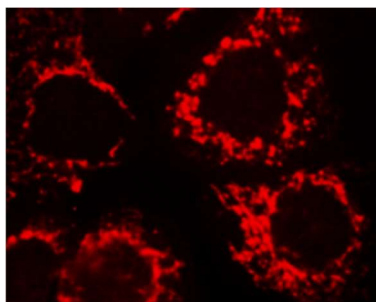

Immunofluorescence staining of formalin-fixed A431 cells showing mitochondrial localization. Detection reagent used: mouse anti-rabbit IgG, PE Conjugated, IRFS4052.

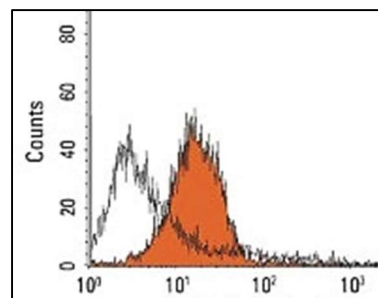

FCM analysis of fixed and permeabilized Raji cells stained with CD45 antibody (IR77-263), followed by PE-conjugated mouse anti-rabbit IgG: IRFS4052. Black line histogram represents the isotype control.

IReal Biotechnology Co., Ltd.

