

Catalog Number: IR53-167

Anti-ZEB1 antibody

Package size: 25, 100 μ l

Store at: -20°C


MW (kDa): 130-200

Overview

Product Name	ZEB1 antibody
Product Number	IR53-167
Gene Description	zinc finger E-box binding homeobox 1
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, mouse, Rat
Recommended Applications Dilutions	Western Blot 1:500 Immunofluorescence 1:200 – 1:300
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	0.68 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

Image

IR53-167 anti-ZEB1 antibody WB image

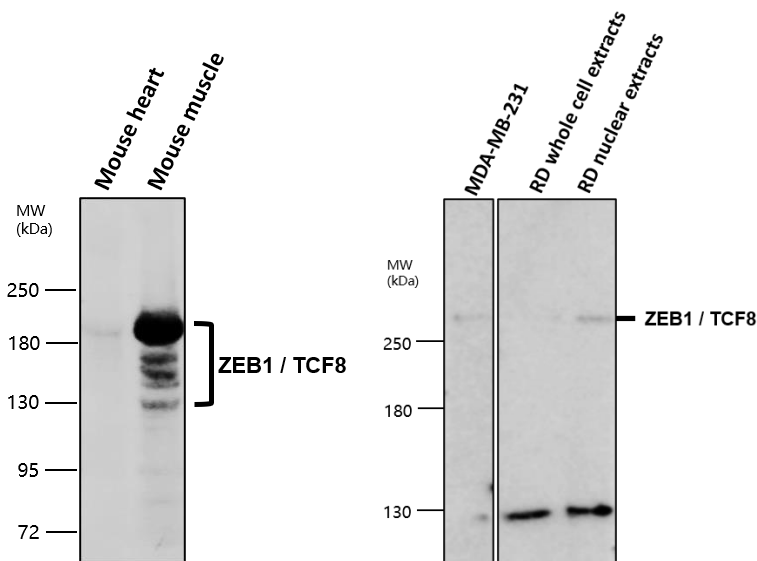

Anti-ZEB1 antibody at 1/500 dilution

Lysates/proteins at 60 μ g per lane
This blot was produced using a 5% SDS-PAGE. The gel was run at 140V for 50 minutes before being transferred onto a Nitrocellulose membrane at 18V for 60 minutes. The membrane was then blocked for an hour before being incubated with IR53-167 overnight at 4°C.