

Catalog Number: IR302-987


Anti-CD81 antibody

Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): 25

Overview

Product Name	CD81 antibody
Product Number	IR302-987
Gene Description	CD81 molecule
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, mouse
Recommended Applications Dilutions	Western Blot 1:1000 – 1:3000 Immunohistochemistry (Paraffin) 1:200 – 1:500 Flow Cytometry 1:100 – 1:200
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	0.75 mg/ml
Purity	Affinity column purified
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR302-987 anti-CD81 antibody WB image

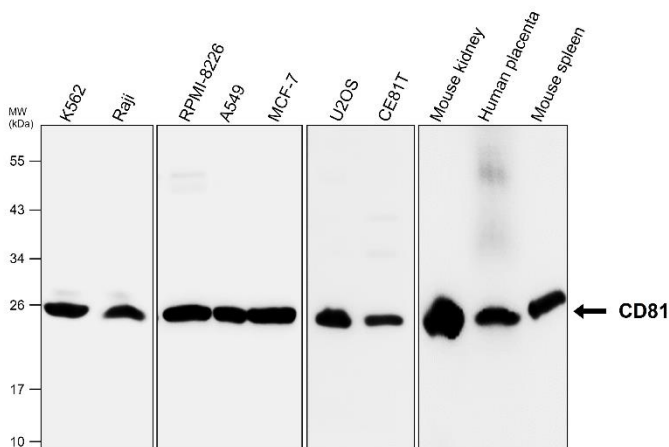

Anti-CD81 antibody at 1/2000 dilution

Lysates at 30 μ g per lane
This blot was produced using a 12% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IR302-987 overnight at 4°C.