

Catalog Number: IR213-777


Anti Estrogen Receptor alpha antibody

Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): 60-66


Overview

Product Name	Estrogen Receptor alpha antibody
Product Number	IR213-777
Gene Description	estrogen receptor 1
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, Mouse, Rat
Recommended Applications Dilutions	Western Blot 1:500 – 1:1000 Immunofluorescence 1:100 – 1:200 Immunohistochemistry (Paraffin) 1:200 – 1:250
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.38 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR213-777 anti- Estrogen Receptor alpha antibody WB image

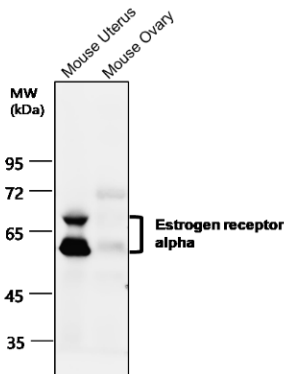

IR213-777 anti- Estrogen Receptor alpha antibody IHC image

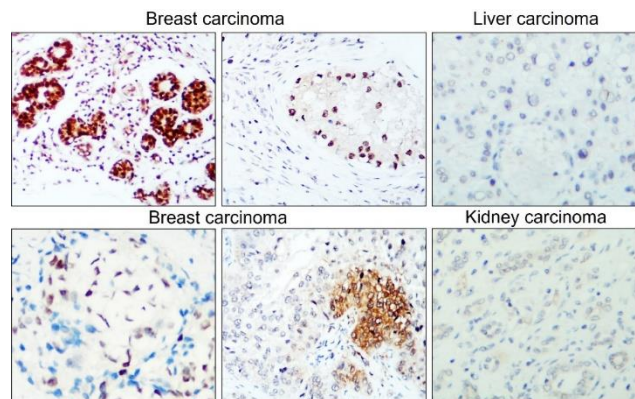

Estrogen Receptor alpha antibody at 1/800 dilution

Lysates at 60 μ g per lane, This blot was produced using a 10% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IR213-777 overnight at 4°C.

Immunohistochemical analysis of paraffin embedded Human cancer tissue labeling Estrogen Receptor alpha with IR213-777 at 1/200.

IReal Biotechnology Co., Ltd.