

WB

IF

IHC

Catalog Number: IRM046


Anti-TSG101, monoclonal antibody

Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): 50


Overview

Product Name	TSG101, monoclonal antibody
Product Number	IRM046
Clonality	Monoclonal
Host	Mouse (mouse IgG1)
Species Reactivity	Human, Mouse, Rat
Recommended Applications Dilutions	Western Blot 1:500 – 1:1000 Immunofluorescence 1:100 – 1:200 Immunohistochemistry (Paraffin) 1:50 – 1:100
Storage Buffer	Phosphate-buffered solution, pH 7.2, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	0.5 mg/ml
Purity	Affinity column purified
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

**IRM046 TSG101, monoclonal antibody
WB image**

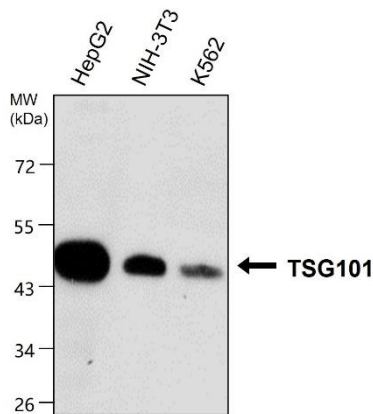

**IRM046 TSG101, monoclonal antibody
IHC image**

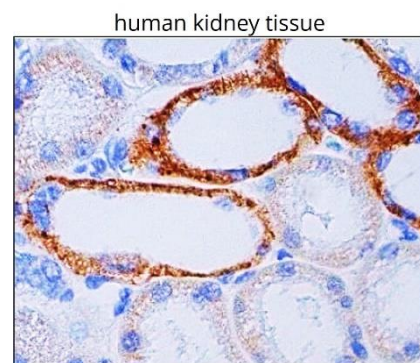

TSG101, monoclonal antibody at 1/800 dilution.
Lysates/ at 35 μ g per lane. Nitrocellulose membrane was then blocked for an hour before being incubated with IRM046 overnight at 4°C.

Immunoperoxidase staining of formalin fixed, paraffin-embedded human kidney tissue.

