

Catalog Number: IRM048


Anti-CDC42 monoclonal antibody

Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): 21


Overview

Product Name	CDC42, monoclonal antibody
Product Number	IRM048
Clonality	Monoclonal
Host	Mouse (mouse IgG1)
Species Reactivity	Human, Mouse, Rat, Drosophila
Recommended Applications Dilutions	Western Blot 1:500 – 1:1000 Immunofluorescence 1:100 – 1:200 Immunohistochemistry (Paraffin) 1:50 – 1:100
Storage Buffer	Phosphate-buffered solution, pH 7.2, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	0.5 mg/ml
Purity	Affinity column purified
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IRM048, CDC42, monoclonal antibody WB image

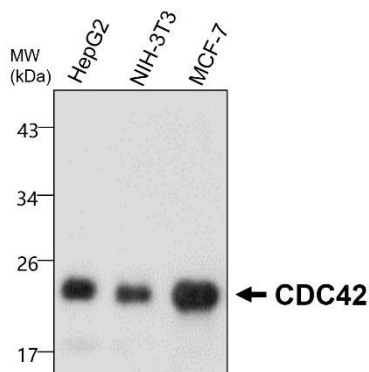

IRM048, CDC42, monoclonal antibody ICC/IF image

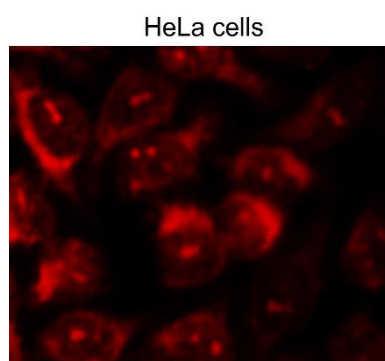

CDC42, monoclonal antibody at 1/500 dilution. Lysates/ at 35 μ g per lane. NC blot was then blocked for an hour before being incubated with IRM048 overnight at 4°C.

Immunofluorescence of HeLa cells.
CDC42, monoclonal antibody (IRM048): 1-200

