

Catalog Number: IR40-117


Anti Phospho-Tau (Thr231) antibody

Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): 52-43


Overview

Product Name	Phospho-Tau (Thr231) antibody
Product Number	IR40-117
Gene Description	microtubule-associated protein tau
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, mouse, rat
Recommended Applications Dilutions	Western Blot 1:1000 – 1:2000 Immunofluorescence 1:200 – 1:400
Storage Buffer	100mM Tris Glycine, 30% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	0.3 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR40-117 anti-Phospho-Tau (Thr231) antibody WB image

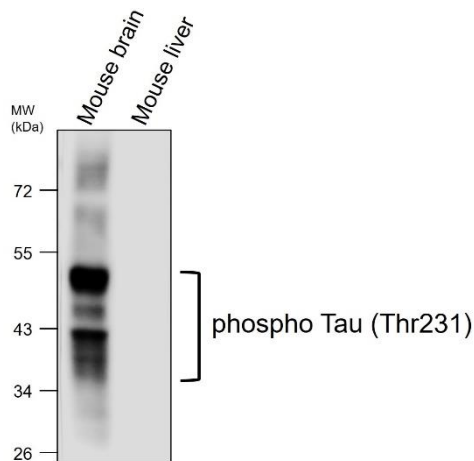

IR40-117 anti-Phospho-Tau (Thr231) antibody WB image

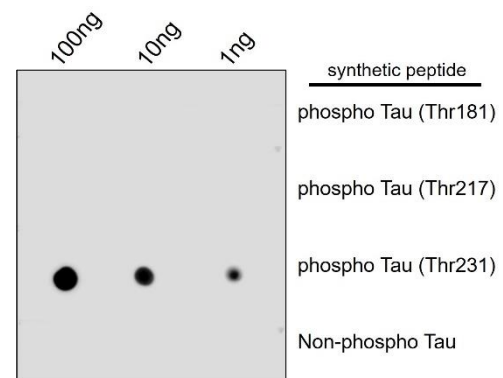

Anti-Phospho-Tau (Thr231) antibody at 1/1000 dilution

Lysates/proteins at 20 μ g per lane
This blot was produced using a 10% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IR40-117 overnight at 4°C.

Dot blot for specificity and sensitivity of **Phospho-Tau (Thr231)** Antibody (IR40-117)

IReal Biotechnology Co., Ltd.