

Catalog Number: IR280-2


## Anti Phospho-EGFR (Tyr1068) antibody

Package size: 25, 100  $\mu$ l  
Store at: -20°C  
MW (kDa): 175


### Overview

Product Name	Phospho-EGFR (Tyr1068) antibody
Product Number	IR280-2
Gene Description	epidermal growth factor receptor
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, mouse, rat
Recommended Applications Dilutions	Western Blot 1:1000 – 1:2000 Immunofluorescence 1:200 – 1:400
Storage Buffer	100mM Tris Glycine, 30% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	0.4 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR280-2 anti-Phospho-EGFR (Tyr1068) antibody WB image

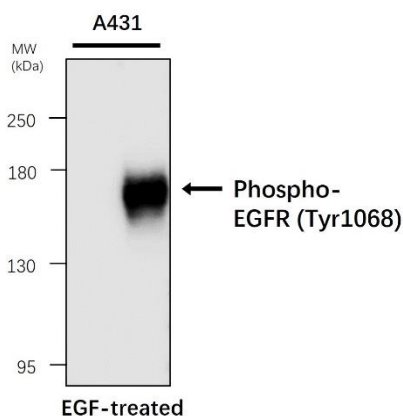

IR280-2 anti-Phospho-EGFR (Tyr1068) antibody WB image

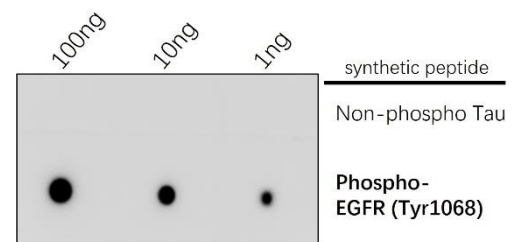

### Anti-Phospho-EGFR (Tyr1068) antibody at 1/1000 dilution

Lysates/proteins at 20  $\mu$ g per lane  
This blot was produced using a 5% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IR280-2 overnight at 4°C.

Dot blot for specificity and sensitivity of Phospho-EGFR (Tyr1068) Antibody (IR280-2)

IReal Biotechnology Co., Ltd.