

Catalog Number: IR283-P1

## Anti Nectin-4 / PVRL4 antibody



Package size: 25, 100  $\mu$ l  
Store at: -20°C  
MW (kDa): 40-50, 55-70

### Overview

Product Name	Nectin-4 / PVRL4 antibody
Product Number	IR283-P1
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human
Recommended Applications Dilutions	Western Blot 1-1000 Immunohistochemistry (Paraffin) 1:100 – 1:200
Storage Buffer	100mM Tris Glycine, 30% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1 mg/ml
Purity	Affinity column purified
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

#### IR283-P1 anti Nectin-4 / PVRL4 antibody IHC image

#### IR283-P1 anti Nectin-4 / PVRL4 antibody WB image

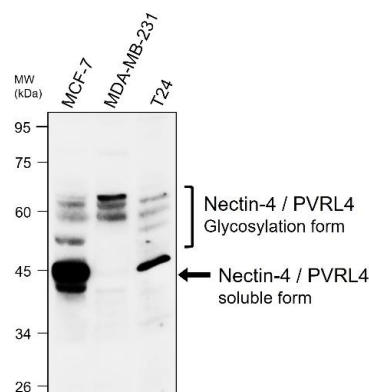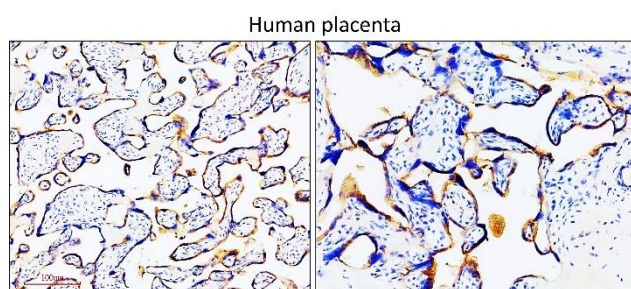

Immunohistochemical analysis of paraffin embedded Human placenta tissue labeling **Nectin-4 / PVRL4** with IR283-P1 at 1/200.

#### Anti Nectin-4 / PVRL4 antibody - at 1/1000 dilution

Lysates/proteins at 25  $\mu$ g per lane  
This blot was produced using a 10% SDS-PAGE. The gel was run at 140V for 50 minutes before being transferred onto a Nitrocellulose membrane at 18V for 60 minutes. The membrane was then blocked for an hour before being incubated with IR283-P1 overnight at 4°C.

